


Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD
 3/4 - 8 Inglewood Place,
 Baulkham Hills NSW 2153 Australia

Test Report Number: 2109108

Testing Requested By: Peter Lee Client's Order Number: Not Supplied

Date Samples Received: 22/09/2021 Date Testing Completed: 8/10/2021

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

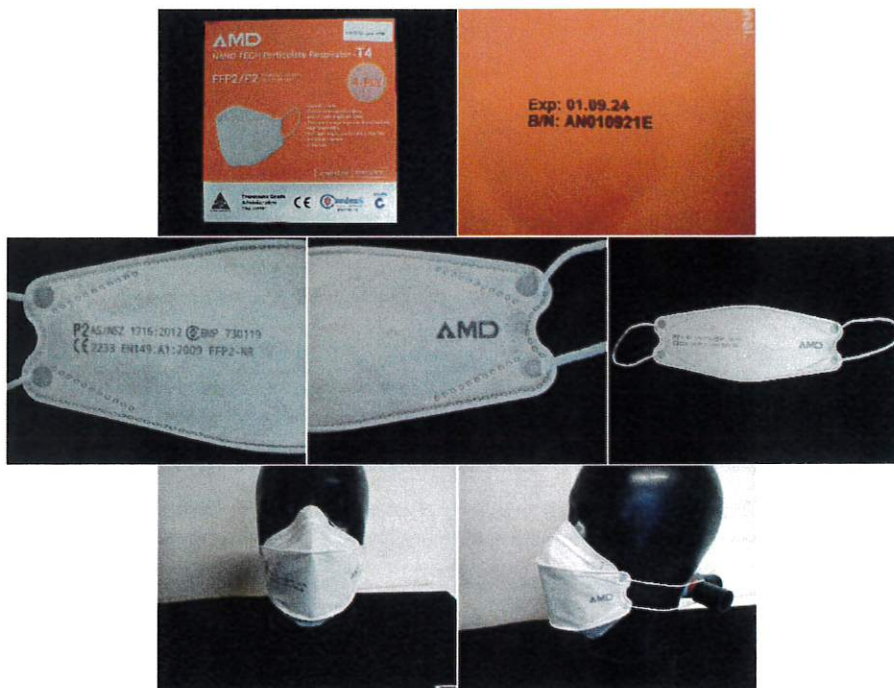

Testing Requested: Class P2
 AS/NZS 1716:2012 – Respiratory Protective Devices, Section 2 – Design and Construction of Assembled Respirators, Clauses 2.1.1, 2.1.2, 2.1.3, 2.1.4, 2.2.1 and 2.2.2 (Appendix D – Total Inward Leakage of Assembled Respirators – Quantitative Sodium Chloride Test)
 AS/NZS 1716:2012 – Respiratory Protective Devices, Section 3 – Facepieces Head Coverings and Harnesses
 AS/NZS 1716:2012 – Respiratory Protective Devices, Section 4 – Particulate Filter Respirators



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:

Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure

R.A. Vickery



Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Pre-Conditioning Summary

Pre-conditioning as per clause 4.3.1 whereby whole filter assembly has been tested and Clause 1.6 has been applied.

Summary

| Clause | Page Number | Compliance/Non-Compliance |
|----------------------|-------------|---------------------------|
| 2.1.1 | 3 | Compliance |
| 2.1.2 | 3 | Compliance |
| 2.1.3 | 4 | Compliance |
| 2.1.4 | 4 | Assessed by BSI |
| 2.2.1 | 5 | Compliance |
| 2.2.2 | 6 | Compliance |
| 3.1.1 | 7 and 8 | Compliance |
| 3.1.4 | 9 | Compliance |
| 3.2.1 | 9 | Compliance |
| 3.2.5 | 10 | Compliance |
| 3.2.6 | 11 | Compliance |
| 4.1 | 11 | Assessed by BSI |
| 4.2.1 | 12 | Compliance |
| 4.2.3 | 12 | Compliance |
| 4.3.3 | 13 | Compliance |
| 4.3.4 | 13 | Compliance |
| 4.3.5 | 14 | Compliance |
| 4.3.6 | 14 | N/A |
| Annex A – Appendix D | 15 | Compliance |

Test Uncertainty:

These uncertainty values are based on a standard uncertainty multiplied by a coverage factor k=2, which provides for a confidence level of approximately 95% - Uncertainty of measurement has not been taken into account when presenting the test result.

NOTE: Samples supplied comply with P2 requirements of AS/NZS 1716:2012 for the clauses tested.



A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure

AUTHORISED SIGNATORY:



R. A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD
 Test Report Number: 2109108
 Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 2 – Design and Construction of Assembled Respirators, Clauses 2.1.1, 2.1.2, 2.1.3, 2.1.4, 2.2.1 and 2.2.2

| Requirement | Result | Compliance / Non-Compliance |
|--|-----------|-----------------------------|
| 2.1 General Requirements | | |
| 2.1.1 Assembled respirators | | |
| Assembled respirators shall be made up of components which have been tested together as a system. | - | N/A |
| The apparatus shall be constructed from durable components | Durable | Compliance |
| The vital parts of the apparatus shall be protected so as to prevent damage and excessive wear during normal use. | Protected | Compliance |
| All parts shall be finished smooth and free from sharp edges and from irregularities that could be a potential hazard or cause discomfort to the wearer. | Smooth | Compliance |

| Requirement | Result | Compliance / Non-Compliance |
|--|-----------------------------------|-----------------------------|
| 2.1.2 Materials | | |
| Respirators should be made of materials able to withstand storage and usage in environments that are likely to be encountered. | Complies | Compliance |
| Materials which may come in contact with the skin should be: | | |
| - Non-staining | Non-staining | Compliance |
| - Soft | Soft | Compliance |
| - Pliable | Pliable | Compliance |
| - Not likely to cause skin irritation | Not likely | Compliance |
| - Shall not taste or smell offensive | Does not taste or smell offensive | Compliance |
| Material from the filtering medium released by air flowing through the filter shall not constitute a hazard or nuisance to the wearer. | Complies | Compliance |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards. The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials. It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so. This document shall not be reproduced except in full

Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

AUTHORISED SIGNATORY:


R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 2 – Design and Construction of Assembled Respirators, Clauses 2.1.1, 2.1.2, 2.1.3, 2.1.4, 2.2.1 and 2.2.2

| Requirement | Result | Compliance / Non-Compliance |
|---|------------|-----------------------------|
| 2.1.3 Filters | | |
| Where appropriate, filters shall be readily replaceable without requiring the use of special tools | Single use | Compliance |
| Shall be designed or marked to prevent incorrect assembly. | One piece | Compliance |
| When the filter is designed to be used with a twin filter facepiece only it shall not be possible to connect the filter to a single filter facepiece unless, by doing so, the respirator assembly would also comply with the requirements for a single filter respirator. | - | N/A |
| The particulate filter of a combined gas and particulate respirator shall be on the influent side of the gas filter. | - | N/A |
| The mass of the replacement filter (or filters) shall not exceed - | - | N/A |
| (a) 300g when it is to be directly connected to a half facepiece; and | - | N/A |
| (b) 500g when it is to be directly connected to a full facepiece. | - | N/A |

| Requirement | Result | Compliance / Non-Compliance |
|---|--------|-----------------------------|
| 2.1.4 Shelf Life | | |
| Each component part of the respirator should have a nominal shelf life of at least five years when properly stored, unless otherwise specified by the manufacturer. | - | Assessed by BSI |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards. The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials. It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so. This document shall not be reproduced except in full

Legend:

NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

AUTHORISED SIGNATORY:


R.A. Vickery



Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 2 – Design and Construction of Assembled Respirators, Clauses 2.1.1, 2.1.2, 2.1.3, 2.1.4, 2.2.1 and 2.2.2

| Requirement | Result | Compliance / Non-Compliance |
|--|--|-----------------------------|
| 2.2 Facial Fit | | |
| 2.2.1 General | | |
| Facial fit of complete respirators shall be tested by determining the total inward leakage of the respirator by a test aerosol of sodium chloride according to the method described in Appendix D. | Complies Size Small A Size Regular B Size Large C | Compliance |
| Personnel shall be selected in accordance with Appendix B. Where the fit of a respirator to a particular physiognomy, whether it be characterized by size or specific facial features, is Assessed by BSI, the testing laboratory, may select test personnel who conform to that physiognomy. Where testing has been carried out in this manner, the manufacturer shall label the respirator accordingly. See Clause 12.1.2.2. | Personnel selected for physiology A, B, C | Compliance |
| When carrying out the test procedure, none of the wearer shall experience any undue discomfort on account of: | | |
| - The fit | Complies | Compliance |
| - Air delivery | Complies | Compliance |
| - Any other feature of the respirator | Complies | Compliance |
| Alternative methods of determining total inward leakage, e.g. using sulfur hexafluoride, may be accepted by the test authority where a correlation with these methods and test criteria has been shown. | - | N/A |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:

Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure



R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD
 Test Report Number: 2109108
 Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 2 – Design and Construction of Assembled Respirators, Clauses 2.1.1, 2.1.2, 2.1.3, 2.1.4, 2.2.1 and 2.2.2

| Requirement | Result | Compliance / Non-Compliance |
|---|---------------------------|-----------------------------|
| 2.2.2 Assessment | | |
| Total inward leakage (TIL) of assembled respirators when tested in accordance with Appendix D shall be assessed for compliance with Table 2.1 | Complies Refer Annex A | Compliance |
| Half facepiece respirators shall be tested with the highest efficiency filters. | - | N/A |
| Full facepiece respirators of the non-powered type shall be tested with either P3 filter or a simulated filter blank. | - | N/A |
| If any test subject records a mean inward leakage for any one test in excess of that specified in Column 2 of Table 2.1, or if the mean of the results of an individual test subject exceeds that specified in Column 1 of Table 2.1, the respirator shall be deemed not to comply with this intended class, or have failed | - | N/A |
| Where the respirator is supplied in more than one size, the test subject shall be supplied with the appropriate size. Any test subject whose maximum inward leakage exceeded that specified in Table 2.1 may participate in further test using an alternative size of the same respirator. | - | N/A |

Table 2.1
 Maximum Total Inward Leakage (TIL) per Test Subject

| Respirator | Percent total inward leakage | |
|--------------------------|---|---|
| | Mean result of test subject not to exceed | No individual exercise result to exceed |
| Non-powered | | |
| - Half face piece | | |
| Class P1 filters | 22.0 | 22.0 |
| Class P2 filters | 8.0 | 8.0 |
| - Full facepiece | | |
| Class P3 | 0.05 | 0.05 |
| Powered | | |
| Class PAPR P1 Filters | 5.0 | 5.0 |
| Class PAPR P2 Filters | 1.0 | 1.0 |
| Class PAPR P3 Filters | 0.05 | 0.05 |
| Air-supplied | | |
| Continuous flow | 0.02 | 0.05 |
| Positive pressure demand | 0.02 | 0.05 |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards. The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials. It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so. This document shall not be reproduced except in full

AUTHORISED SIGNATORY:



R.A. Vickery



Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD


Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 3 – Facepieces Head Coverings and Harnesses

| Requirement | Result | Compliance / Non-Compliance |
|---|---------------------|-----------------------------|
| 3.1 DESIGN REQUIREMENTS | | |
| 3.1.1 GENERAL | | |
| Each facepiece or head covering shall comply with the following: | | |
| (a) Be designed to fit a wide range of facial contours and head sizes of the workplace population | Standard Size | Compliance |
| (b) Be supported on the head or shoulders by suitable means so that the device remains in position during normal work practices, especially when the wearer bends forward from the waist | Remains in position | Compliance |
| (c) Permit the component parts likely to require service to be readily detached for maintenance and cleaning, but be secure against accidental disconnection. | - | N/A |
| (d) Where the head covering has been designated by the manufacturer as being suitable for abrasive blasting, the construction shall provide physical protection to the wearer's head, shoulders and upper part of the chest against rebounding abrasives. | - | N/A |
| Helmets fitted with a visor which is supported in a hinged frame shall have a means for securely fastening the frame in its closed position so that it cannot be opened inadvertently. | - | N/A |

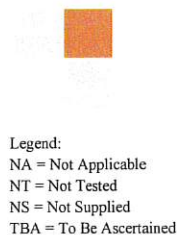

Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:



R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 3 – Facepieces Head Coverings and Harnesses

| Requirement | Result | Compliance / Non-Compliance |
|--|----------|-----------------------------|
| (e) Where an inner bib is proved as part of the head covering, it shall be of a material which will prevent or restrict the flow of air through it | - | N/A |
| The bib may have a drawstring or elasticized neck band or cuff to draw the bib closely around the wearers neck | - | N/A |
| The hood or bib shall be readily removable for cleaning or replacement. | - | N/A |
| The design of the facepiece or head covering should cause the least possible interference with speech and vision. | Complies | Compliance |
| The full facepiece should be designed to minimize misting of the face mask, e.g. by provision of orinasal inserts or nose cups. It should also permit the use of special spectacles designed for use without temple pieces, so that air tightness is not affected. | - | N/A |
| Full facepieces should incorporate facilities for speech transmission. The components of any electrically operated speech transmission device shall be 'intrinsically safe' or 'flameproof' (see Clause 2.1.7) if they are to be used in flammable atmospheres. | - | N/A |
| Where the use of the apparatus is intended solely for escape, a mouthpiece an nose clip may be incorporated in place of a facepiece | - | N/A |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:

Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure



R.A. Vickery



Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 3 – Facepieces Head Coverings and Harnesses

| Requirement | Result | Compliance / Non-Compliance |
|---|----------|-----------------------------|
| 3.1.4 Nose Clip | | |
| The nose clip shall be designed so as to effect the maximum possible security against displacement or slipping, for example, if the wearer receives a chance blow, or stumbles, or their nose becomes wet with perspiration | Complies | Compliance |
| Suitable means shall be provided for attaching the nose clip to the mouthpiece | Suitable | Compliance |
| The nose clip shall be so positioned that, when inserting the mouthpiece, the user is made aware of the need to apply the nose clip | Complies | Compliance |

| Requirement | Result | Compliance / Non-Compliance |
|--|----------|-----------------------------|
| 3.2 Performance Requirements | | |
| 3.2.1 Facial fit | | |
| In combination with other components, e.g. filters and air supply, the assembled respirator shall provide adequate protection either by means of a facial seal or by the provision of positive pressure in the space enclosed by the respirator, or by both, to minimize the entry of ambient atmosphere | Complies | Compliance |
| Facial fit of complete respirators shall be tested by determining the total inward leakage in accordance with Clause 2.2 | Complies | Compliance |



A2LA Accredited Laboratory
 Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards. The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials. It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so. This document shall not be reproduced except in full

Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

AUTHORISED SIGNATORY:



R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108


Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 3 – Facepieces Head Coverings and Harnesses

| Requirement | Result | Compliance / Non-Compliance |
|--|--------|-----------------------------|
| 3.2.4.2 Exhalation valve leakage | | |
| The exhalation valve or fitting surface of the facepiece shall be sealed to a former and, where appropriate, the inlet air opening of the facepiece sealed. | - | N/A |
| When tested with air at a constant suction head of 250Pa, the leakage into the facepiece from the valve or valves and any leakage as appropriate from around the eyepiece or visor shall not exceed 30 mL/min. | - | N/A |
| During the test, both valve(s) and seating shall be free of moisture | - | N/A |

| Requirement | Result in Pa | | | | | Compliance / Non-Compliance |
|---|---------------|-----|----|-----|----|-----------------------------|
| | Sample Number | | | | | |
| 3.2.5 Exhalation resistance – Air filtering respirators | | | | | | |
| When tested in accordance with Appendix G at a continuous flow rate of 85 ± 2 L/min, the exhalation resistance of the entire assembly, measured relative to the static pressure in the facepiece, shall be less than or equal to- | 1 | 2 | 3 | 4 | 5 | |
| (a) For all full facepieces.....200 Pa | - | - | - | - | - | N/A |
| (b) For all half facepieces.....120 Pa | 104 | 107 | 98 | 100 | 97 | Compliance |

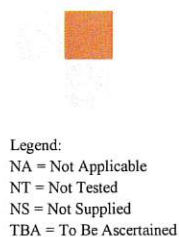

Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards. The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials. It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so. This document shall not be reproduced except in full

AUTHORISED SIGNATORY:



R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 4 – Particulate Filter Respirators

| Requirement | Result | Compliance / Non-Compliance |
|--|----------|-----------------------------|
| 3.2.6 Security of attachments | | |
| Fittings directly attached to the head covering or facepiece including filter receptacles, exhalation valve housing, speech diaphragms and demand valves but excluding straps and buckles shall be tested accordingly for security of attachment as applicable when assembled in accordance with the manufacturer's instructions | - | N/A |
| Each such attachment shall withstand an axial tensile force of 50N for 10s | - | N/A |
| Each strap, buckle and its attachment to a half facepiece shall withstand an axial tensile force of 10 N applied for 10s in the direction of pulling when the facepiece is fitted. | Complies | Compliance |

| Requirement | Result | Compliance / Non-Compliance |
|---|--------------------------------------|---|
| 4.1 Design and Construction | | |
| The design and construction of a particulate filter shall be such that, when combined or incorporate with the appropriate facepiece or head covering it shall provide protection against particulates in accordance with its class. | | |
| The respirator shall be designed so that all the inhaled air passes through the filter or filters. | Complies | Compliance |
| Particulate filter respirators shall comply with Sections 2, 3 and 12, as appropriate. A filter shall comply with Clause 4.3 | Section 2 Section 3 Section 12 | Compliance Compliance Assessed by BSI |



A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:
Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure



R.A. Vickery



CERT #6187.01

Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:
AS/NZS 1716:2012 – Respiratory Protective Devices, Section 4 – Particulate Filter Respirators

| Requirement | Result | Compliance / Non-Compliance |
|---|-------------|-----------------------------|
| 4.2 Classification and Components | | |
| 4.2.1 Classes | | |
| Three classes of particulate filter are differentiated according to filtering efficiency: | | |
| (a) Class P1 – intended for use against mechanically generated particulates of sizes most commonly encountered in industry. | - | N/A |
| (b) Class P2 – intended for use against both mechanically and thermally generated particulates. | P2 complies | Compliance |
| (c) Class P3 – intended for use against all particulates including highly toxic materials | - | N/A |

| Requirement | Result | Compliance / Non-Compliance |
|---|----------------|-----------------------------|
| 4.2.3 Components | | |
| The respirator shall include a particulate filter complying with Clause 4.3 | | |
| It may incorporate - | | |
| (a) A full or half facepiece, a head covering, or a mouthpiece held securely in position by a head harness; | Half facepiece | Compliance |
| (b) An exhalation valve or exhalation valve assembly; | - | N/A |
| (c) An inhalation valve or inhalation valve assembly; | - | N/A |
| (d) One or more filter holders; | - | N/A |
| (e) One or more flexible breathing tubes; | - | N/A |
| (f) A belt or harness to attach the filter or filters to the wearer's body; and | - | N/A |
| (g) A gas and vapour filter complying with Clause 5.4 | - | N/A |


A2LA Accredited Laboratory
Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:
Legend:

NA = Not Applicable

NT = Not Tested

NS = Not Supplied

TBA = To Be Ascertained

Document name: Differential Pressure


R.A. Vickery


Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:


AS/NZS 1716:2012 – Respiratory Protective Devices, Section 4 – Particulate Filter Respirators

| Requirement | Result | Compliance / Non-Compliance |
|---|------------|-----------------------------|
| 4.3.3 Simulated Wear Treatment | | |
| Respirators except PAPR which have no exhalation valve, or where a substantial proportion of the exhaled air passes back through the filter, shall be subject to exhaled air humidity pre-conditioning in accordance with Paragraph E5.6 of Appendix E. | Subjected | Compliance |
| If filters are not tested directly upon completion of such pre-conditioning, they should be stored in a manner so as to retain their humidity prior to testing for filtering efficiency. | Sealed Bag | Compliance |

| Requirement | Result in Pa Sample Number | Compliance / Non-Compliance |
|---|-----------------------------|-----------------------------|
| 4.3.4 Inhalation resistance | | |
| When tested in accordance with Appendix G with a continuous stream of air passing through the assembly at a defined rate, the resistance imposed by the assembly shall not exceed the values given in Table 4.1 | | |
| | 1 2 3 4 5 | |
| at 30 ± 1 L/min | 41 42 37 39 39 | Compliance |
| at 95 ± 2 L/min | 143 144 138 140 139 | Compliance |

 TABLE 4.1
 INHALATION RESISTANCE

| Filter class | Filter assembly only maximum resistance, Pa* | | Assembled respirator maximum resistance, Pa* | |
|--------------|--|----------------|--|----------------|
| | At 30 ±1 L/min | At 95 ±2 L/min | At 30 ±1 L/min | At 95 ±2 L/min |
| P1 | 60 | 210 | 110 | 340 |
| P2 | 70 | 240 | 120 | 370 |
| P3 | 120 | 420 | 170 | 570 |

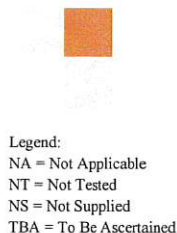* 1 mbar = 100 Pa = 100 mm H₂O


Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:


R. A. Vickery



Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD
Test Report Number: 2109108
Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

Results:

AS/NZS 1716:2012 – Respiratory Protective Devices, Section 4 – Particulate Filter Respirators

| Requirement | Result in % <small>Per TGA 10 December 2020 Evidence requirements for face masks that are medical devices</small> | | Compliance / Non- Compliance |
|---|--|------|------------------------------------|
| 4.3.5 Test of filtering efficiency | | | |
| When sealed to a suitable former, and tested in accordance with Appendix I, non-powered respirator filters shall not show penetration in excess of the following: | | | |
| (a) Class P2.....not more than 6% | Specimen No.1 | 0.26 | Compliance |
| | Specimen No.2 | 0.27 | Compliance |
| | Specimen No.3 | 0.22 | Compliance |
| | Specimen No.4 | 0.24 | Compliance |
| | Specimen No.5 | 0.35 | Compliance |
| | Specimen No.6 | 0.24 | Compliance |
| | Specimen No.7 | 0.26 | Compliance |
| | Specimen No.8 | 0.22 | Compliance |
| | Specimen No.9 | 0.28 | Compliance |
| | Specimen No.10 | 0.24 | Compliance |
| When a single filter of a twin filter respirator is tested separately, the air flow specified for this test shall be halved. If it is possible that the single filter may be used in a single filter respirator, then the full airflow shall be used. | Single Filter | | N/A |

| Requirement | Result | Compliance / Non-Compliance |
|--|--------|--------------------------------|
| 4.3.6 Filters used in series | | |
| Where a separate particulate filter is used in series with any other filter, the particulate penetration shall be tested in the combined configuration in accordance with Clause 4.3.5 | - | N/A |

Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

A2LA Accredited Laboratory
Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing
 The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.
 The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.
 It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.
 This document shall not be reproduced except in full

AUTHORISED SIGNATORY:

R. A. Vickery



Client: ADVANCED MEDICAL DEVICES (AMD) PTY LTD

Test Report Number: 2109108

Sample Description: New T4 printed logo and text, labelled Nano P2 earloop, Client advises Batch number AN010921E, Half facepiece, Tri-fold particulate respirator, Internal noseclip, Internal nose foam padding, Twin earloops, Printed marking, Samples as supplied.

ANNEX A


Results:
 AS/NZS 1716:2012 – Respiratory Protective Devices, Clauses 2.2.1 and 2.2.2 (Appendix D – Total Inward Leakage of Assembled Respirators – Quantitative Sodium Chloride Test)

Total Inward Leakage (TIL) %

| Test Subject (Person) | Facial Dimension Group (As per Appendix B) | Walk | Side to Side | Up/Down | Talk | Walk | Mean |
|-----------------------|--|------|--------------|---------|------|------|------|
| 1 | A | 4.45 | 6.36 | 5.18 | 3.27 | 4.34 | 4.72 |
| 2 | A | 4.21 | 5.16 | 4.02 | 4.00 | 5.23 | 4.52 |
| 3 | B | 5.11 | 5.22 | 6.88 | 4.11 | 6.48 | 5.56 |
| 4 | B | 4.87 | 5.87 | 6.48 | 3.25 | 4.61 | 5.02 |
| 5 | B | 5.02 | 4.12 | 5.21 | 2.01 | 5.78 | 4.43 |
| 6 | B | 4.36 | 3.88 | 5.05 | 2.99 | 5.34 | 4.32 |
| 7 | B | 3.99 | 6.31 | 6.45 | 3.21 | 3.88 | 4.77 |
| 8 | B | 4.11 | 7.55 | 6.22 | 3.77 | 4.25 | 5.18 |
| 9 | C | 4.77 | 3.61 | 6.01 | 2.89 | 4.36 | 4.33 |
| 10 | C | 5.00 | 6.33 | 6.47 | 3.06 | 5.20 | 5.21 |
| Mean | | | | | | | 4.81 |

| Class P1 Filters – Requirement | Result | Compliance / Non-Compliance |
|---|--------|-----------------------------|
| Individual exercise results not to exceed 22% (as per Table 2.1) | - | N/A |
| Mean result of Test Subjects not to exceed 22% (as per Table 2.1) | - | N/A |

| Class P2 Filters – Requirement | Result | Compliance / Non-Compliance |
|--|----------|-----------------------------|
| Individual exercise results not to exceed 8% (as per Table 2.1) | Complies | Compliance |
| Mean result of Test Subjects not to exceed 8% (as per Table 2.1) | Complies | Compliance |



Legend:
 NA = Not Applicable
 NT = Not Tested
 NS = Not Supplied
 TBA = To Be Ascertained

Document name: Differential Pressure

A2LA Accredited Laboratory

Certificate Number 6187.01

Accredited for compliance with ISO/IEC 17025 - Testing

The results of the tests, calibrations and/or measurements in this document are traceable to Australian National Standards.

The results are applicable to the sample tested and may not apply to other batches of the same material or similar materials.

It is a condition of the provision of these test results that you do not use the name of VicLab in connection with the promotion or sale of your products, unless VicLab has given express written authority to do so.

This document shall not be reproduced except in full

AUTHORISED SIGNATORY:



R. A. Vickery

